



PURE**RELIEF LUXE**

MICROMINK HEATING PAD



USER MANUAL MODEL: PELUXWAV-S

IMPORTANT SAFETY INSTRUCTIONS DO NOT DESTROY

DANGER: TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE **FOLLOWING INSTRUCTIONS:**

- READ ALL INSTRUCTIONS.
- DO NOT USE WHILE SLEEPING.
- DO NOT USE ON AN INFANT OR ON AN ANIMAL.
- THIS PRODUCT IS NOT TO BE USED ON OR BY AN INVALID, A SLEEPING OR UNCONSCIOUS PERSON, PERSON WITH DIABETES, OR PERSON WITH POOR BLOOD CIRCULATION UNLESS CAREFULLY ATTENDED. DO NOT USE ON AREAS OF THE SKIN THAT ARE
- INSENSITIVE TO TOUCH OR HEAT.
- BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING. CHECK SKIN UNDER PRODUCT FREQUENTLY.
- DO NOT USE IN AN OXYGEN-FILLED ATMOSPHERE.
- DO NOT USE PINS OR OTHER METALLIC OBJECTS TO FASTEN THIS PRODUCT INTO PLACE.
- DO NOT SIT ON OR AGAINST THE PRODUCT OR CRUSH IT. AVOID SHARP FOLDS IN THE FABRIC.
- 10. DO NOT USE THE CORD AS A HANDLE.
- 11. CHILDREN SHOULD NOT USE THIS UNIT WITHOUT ADULT SUPERVISION.
- 12. PLACE THE UNIT ON TOP OF AND NOT UNDER THE PART OF THE BODY THAT REQUIRES HEAT THERAPY.

- 13. ONLY USE THIS PRODUCT ON A 120V AC CIRCUIT.
- 14. UNPLUG WHEN NOT IN USE.
- 15. DO NOT USE WITH GENERATORS, POWER CONVERTERS, OR INVERTERS.16. DO NOT TAMPER WITH THIS PRODUCT IN ANY WAY.
- THERE ARE NO USER-SERVICEABLE PARTS. IF FOR ANY REASON THIS PRODUCT DOES NOT FUNCTION PROPERLY, CONTACT PURE ENRICHMENT AT HELP@PUREENRICHMENT.COM

 17. DO NOT USE THIS PRODUCT WITH LINIMENT, SALVE
- OR OINTMENT PREPARATIONS AS SKIN BURNS CAN RESULT.
- 18. NEVER PULL THIS PAD BY THE SUPPLY CORD.
- 19. DO NOT BEND OR PINCH THE CORD.
- 20. DO NOT USE THE PRODUCT IF THE CORD IS DAMAGED.
- 21. NEVER LEAVE THE PRODUCT UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
- 22. LOOP CORD LOOSELY WHEN STORING, TIGHT WRAPPING CAN DAMAGE CORD AND INTERNAL PARTS.
- 23. THIS PRODUCT HAS A POLARIZED PLUG (ONE BLADE IS WIDER THAN THE OTHER). TO REDUCE THE RISK OF ELECTRIC SHOCK, THIS PLUG IS INTENDED TO FIT IN A POLARIZED OUTLET ONLY ONE WAY. IF THE PLUG DOES NOT FIT FULLY IN THE OUTLET, REVERSE THE PLUG. IF IT STILL DOES NOT FIT, CONTACT A QUALIFIED ELECTRICIAN TO INSTALL THE PROPER OUTLET.

DO NOT ATTEMPT TO MODIFY THE PLUG OR DISABLE THIS SAFETY FEATURE IN ANY WAY.


24. SAVE THESE INSTRUCTIONS.

INTRODUCTION

Thank you for purchasing PureRelief Luxe. We designed this plush heating pad to provide temporary soothing heat therapy to large muscle groups in your back, legs, arms, and more.

USING PURERELIEF LUXE

Attach the Pad Connector to the Control Connector. Push both parts together securely to ensure a tight fit.

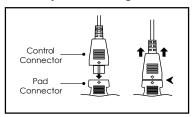

NOTE: The Control Connector must be connected to the Pad Connector before plugging the unit into an electrical outlet.

- After connecting the Pad Connector to the Control Connector, plug the unit into an electrical outlet.
- To turn PureRelief Luxe on, press the power button. Use the Heat ▲ and Heat ▼ buttons to select the desired heat setting:

1=Warm 4=Medium 2=Low 5=Medium High

3=Medium Low 6=High

To turn PureRelief Luxe off, press the power button one more time.

NOTE: If the unit does not turn on when you press the power button, make sure the Control Connector is securely plugged into the Pad Connector

AUTO-OFF FEATURE

To conserve energy and prevent skin burns caused by excessive use, PureRelief Luxe includes an auto-off feature. The unit will turn off automatically after two hours of continuous use. You can turn this feature off by pressing the Auto Off button until the Stay On indicator light turns red. The unit will not turn off automatically after two hours of continuous use when the Stay On light is red.

NOTE: After a power outage, PureRelief Luxe automatically resets the timer on the auto-off feature.

WASHING INSTRUCTIONS

DO NOT DRY CLEAN! DO NOT USE BLEACH! DO NOT USE WRINGER!

- Disconnect PureRelief Luxe from the electrical outlet.
- Detach the Pad Connector from the Control Connector.
- Machine wash PureRelief Luxe in cold water on the gentle cycle, or hand wash. Line dry or tumble dry on high heat for five minutes. If using an automatic dryer, remove while damp and drape to finish drying.
- 4. Do not iron the heating pad.
- Do not plug in PureRelief Luxe until the Pad Connector is completely dry.

PRODUCT SPECIFICATIONS

Pad Size:	12" x 24"
Power Supply:	AC: 120V, 60Hz, 180W
Pad Material:	Luxury Micromink (Top); Polyester (Bottom)
Auto Shut-Off Timer:	120 minutes
6 Heat Settings:	Warm–temperature range (~105°F) Low–temperature range (~112°F) Medium Low–temperature range (~119°F) Medium–temperature range (~126°F) Medium High–temperature range (~133°F) High–temperature range (~140°F)

WARRANTY

PureRelief Luxe comes with an industry-leading 5 Year Warranty that begins on the date of purchase.

The warranty applies to the heating pad and necessary parts and labor relating thereto. The warranty does not apply to damage resulting from failure to follow the operating instructions, accidents, abuse, alterations, or disassembly by unauthorized individuals.

PRODUCT REGISTRATION

To ensure your full warranty and to receive product updates and streamlined customer support, remember to register your product at:

www.pureenrichment.com

ARE YOU 100% SATISFIED?

If you have product questions or concerns don't hesitate to contact us directly at:

help@pureenrichment.com

Pure Enrichment has an award-winning warranty, exchange, and customer service program that guarantees hassle-free solutions within 24 hours to any issue you might have!

FCC NOTICE

Model: PureRelief Luxe Micromink Heating Pad (PELUXWAV-S)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Pure Enrichment 14505 Astronautics Ln Huntington Beach, CA 92647 www.pureenrichment.com

The product is protected against unauthorized use. DO NOT modify it without the authorization of the manufacturer, or the user's authority to operate could be voided.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.



Huntington Beach, CA 92647 www.pureenrichment.com