

ITEM # OBT24

OVER BED TABLE 24" x 16"

ASSEMBLY & OPERATING INSTRUCTIONS



READ ALL INSTRUCTIONS AND WARNINGS BEFORE USING THIS PRODUCT

This manual provides important information on proper operation and maintenance. Every effort has been made to ensure the accuracy of this manual. We reserve the right to change this product at any time without prior notice. IF THERE IS ANY QUESTION ABOUT THIS PRODUCT BEING SAFE OR UNSAFE, DO NOT USE THIS PRODUCT!

QUESTIONS? PROBLEMS? CALL CUSTOMER SERVICE AT 1-636-532-9888

If you need instructions on how to use this item properly, experience a problem, or need parts, call our customer service department Monday-Friday, 8 AM - 4 PM Central Time. A copy of the sales receipt is required.


WARNING: DO NOT OVERLOAD THE MAXIMUM WEIGHT CAPACITY

This unit can support a maximum weight capacity of 22 lbs.

ASSEMBLY


STEP 1

Use hex wrench to attach casters to H-shaped sub-frame (FIG 1)

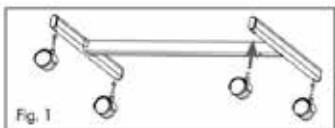

STEP 2

Attach table frame to H-shaped sub-frame by using the two washers plus the screws. Attach to the central pole. (FIG 2)

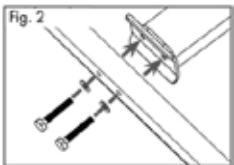


STEP 3

Press in the metal pin while removing the plastic cap. (FIG 3)



STEP 4

Push table frame through the eyelets which are fastened to the tabletop. At the same time, connect the ring which is fastened to the tabletop for incline setting to the central pole. (FIG 4)

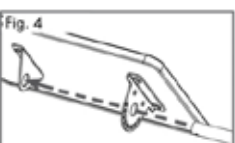

READ ALL INSTRUCTIONS BEFORE ASSEMBLY

DO NOT EXCEED 22 LB WEIGHT CAPACITY

MAKE SURE ALL CASTERS ARE LOCKED BEFORE USE.

CONTENTS

4 PIECES
LOCKING CASTERS



1 PIECE
SUB-FRAME



2 PIECES
SCREWS



1 PIECE
TABLE TOP



2 PIECES
WASHERS



1 PIECE
STAND PIPE



1 PIECE
HEX WRENCH

WARNING!

THIS IS NOT A TOY. ONLY USE AS INTENDED BY ADULTS

DO NOT EXCEED 22 LB WEIGHT CAPACITY

MAKE SURE ALL CASTERS ARE LOCKED BEFORE USE.

MAKE SURE ALL SCREWS ARE TIGHT BEFORE USE.

TABLE CAN TIP OVER WHEN CHILDREN ATTEMPT TO PULL

THEMSELVES UP ON IT, OR SIT ON THE TABLE

DO NOT SIT OR STAND ON TABLE.


STEP 5

Push the table frame through both eyelets, and replace the plastic cap on the end of the table frame. (FIG 5)




ADJUSTING THE TABLE HEIGHT

Loosen the screws with the plastic grip and set the table to the desired height. Then retighten the screws. (FIG 6)

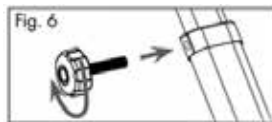

ADJUSTING THE TABLE TILT ANGLE

Pull the locking pin out and adjust the table tilt angle. Return the locking pin to its original position. (FIG 7)

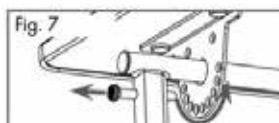

CARE, STORAGE AND MAINTENANCE

Only use warm water or mild cleaning products and soft sponge to clean the surfaces. Always check the screw connections periodically.